

DilaSoft I version 2.75


Windows Software for Theta Dilatometers; Rotating, Parallel Plate, Bending Beam Viscometers and DTA/TGA

Theta's new Windows software features multiple larger and improved **interactive run-time graphics displays** in a **multitasking** environment allowing experiments to run while simultaneously analyzing previous experiment data.

The run-time display shows all important information as the experiment progresses. The displays include: real-time bar graphs of setpoint and specimen temperatures, numeric and strip chart recorder of Measured Values (MV)* and temperature, temperature segment profile and much more. While the experiments are running, interactive graphs can be displayed showing MV* vs. temperature and/or MV* vs. time vs. temperature. These graphs can be moved, sized and scaled while the experiment continues to run, giving the operator all necessary information about the experiment as it is taking place.


The **Quick Report Generator** automatically gathers all the data about the experiment including the materials, temperature program, tables, dates, graphs and data and formats them into user selectable reports that can be saved to disk files for later use or printed immediately on the default Windows printer.

* Measured Value defined by instrument type.
 Dilatometer = Expansion. Viscometer = Viscosity. DTA = Temperature difference. TGA = mass. Etc.

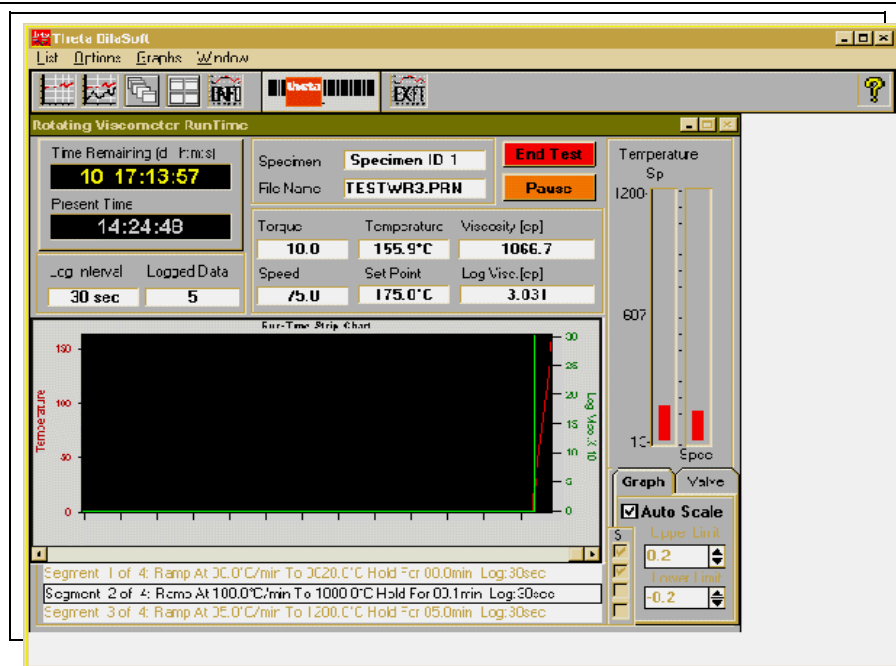

Data analysis remains unsurpassed in the industry offering the widest array of tools, charts, graphs, tables, printouts and reports. Analysis graphs and reports can all be completely customized by the user, printed immediately or saved to disk files for later retrieval or imported into other applications.

Theta Industries has been the industry leader in thermal testing for over 30 years. We remain dedicated to bringing our customers the finest in both software and hardware to create the most advanced thermal testing equipment available.

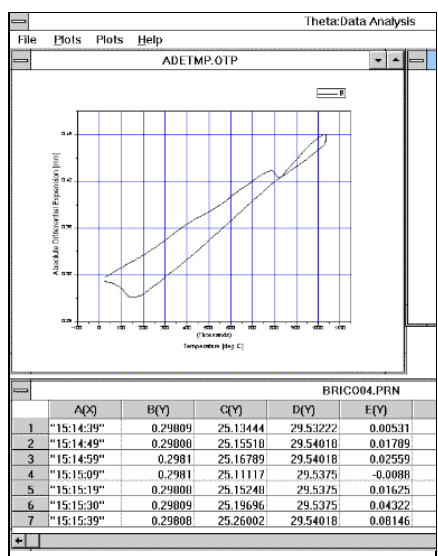

For information about the price and availability of Theta's latest software and equipment offerings, call, write, FAX or E-mail to the address below for a prompt response. Or see our web site for the latest news, information and demos.



DilaSoft I version 2.75

New Features

- ♦ Familiar Windows / mouse interface with larger, easier to read windows complete with full on-line Windows help files, accessible at any time.
- ♦ Channel Mapping allows setting the origin and type of any channel.
- ♦ Correction files automatically correct for any variance in equipment for exact measurements.
- ♦ Standard reference materials for Viscometers and Dilatometers plus user entered reference materials.
- ♦ "Change-Over" height & deflection is a user setting that forces the software to jump to the next temperature program segment when reached, for Parallel Plate and Bending Beam Viscometers..
- ♦ Pre-heat and Delayed Start.
- ♦ Improved temperature profile programming 100 segments – each with its own data storage interval, ramp rate, hold time, digital switch and/or mass flow valve control, head speed (Rotating Viscometer), text alert, and more. The entire profile is displayed graphically on screen and can be printed.
- ♦ Optional use of maximum length for shrinkage calculation, valuable to glass analysis (Dilatometer).
- ♦ Automatic data file naming for fast repetitive experiments.
- ♦ Specimen length auto-length for faster experiment starts without having to measure the specimen (Dilatometer).
- ♦ Optional user selectable reset temperature for safer experiment shutdown or leave the furnace at a high temperature for repeated experiments..
- ♦ Instrument calibration is fully controlled through software. A selectable calibration interval reminds the operator when the equipment requires calibration
- ♦ Includes Setpoint utility diagnostic program..

Quick Reports

Automatic reports:

- ♦ Summary of experiment data
- ♦ Full temperature program
- ♦ Table of results
- ♦ Full data file

Automatic charts:

- ♦ MV* vs. temperature
- ♦ Average alpha* vs. temperature
- ♦ MV* vs. time
- ♦ Temperature vs. time

Run-time Display

- ♦ Audible/visual alerts and messages (the software actually speaks, with optional hardware). These can be changed and set by the user.
- ♦ Rotating Viscometer is now allowed greater than 100% torque spikes without shutting down the experiment.
- ♦ View or edit the entire temperature program while the experiment runs.
- ♦ Experiment parameters can be changed while the experiment runs.
- ♦ Graphical LVDT positioning for added precision (Dilatometer).
- ♦ Simulated strip-chart recorder displays the last 750 measured values and temperature points in a scrollable, time-based window.
- ♦ Live displays of measured and calculated values; the current, prior and next temperature program segments; setpoint and actual specimen temperature. A mini-Status window is available to save screen space.
- ♦ Icon-style toolbar permits easy selection of graphs and overall window organization.
- ♦ Run-time graphs include MV* vs. temperature, and MV* and temperature vs. time. Each graph is an independent window which can be moved, sized, printed or closed while the experiment continues executing in the background.

* Value depends on instrument type.

Data Analysis

- ♦ Experiments can be run in the background while analysis is performed in the foreground.
- ♦ Floating toolbar for graph operations (zoom, point selection, etc.).
- ♦ Multiple windows for graphs and data.
- ♦ Many pre-formatted graph templates supplied for user convenience, each set corresponding to specific instrument type.
- ♦ Automatic TMA and peak detection functions.
- ♦ User selectable reference/offset temperature, decimal points and temperature scale (C, F and K).
- ♦ Adjustable data smoothing.
- ♦ Automatic temperature transition point determination and plotting.
- ♦ Complete *Table of Results* with user selectable temperature range and increment parameters.
- ♦ All reports can be viewed, printed or stored as disk files for later retrieval or importing into other applications.
- ♦ Data files are now registered with the Analyze program so they can be double-clicked to open.

Theta Industries, Inc. is committed to providing our customers the leading edge of technology, service and support with the finest instruments available.

There are many new features and updates to our software.

For the latest update information about DilaSoft I please see our web page at:

www.theta-us.com/software/DilaSoft1.html

Phone, fax or write for further information, pricing or upgrade options. Or see our Web site for more details, the latest news and demos.